FREEQUENT

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Freequent™

Radiofrequency and cooling beauty device

- Increases skin radiance, density
- Boosts collagen naturally
- Brings professional level aesthetic results to your home
- Reduces wrinkles and fine lines

User Guide

Important Safety Warnings!

-  Children should be supervised at all times to ensure they do not play with the device.
-  Please read all warnings and safety precautions before use, and strictly follow the guidelines when using Freequent™.
-  Keep away from the children. Close supervision is necessary when this product is used by, on, or near children or physically challenged persons.
-  Not servicing/maintenance while the equipment is in use.
-  The appliance is only to be used with the unit provided.
-  Children shall not play with the appliance.
-  Cleaning and user maintenance shall not be made by children without supervision.
-  **Warning !**
Using the Qlarite "Freequent™" device without the Freequent Radiance Enhancing RF Gel or using any other gel or material, will adversely affect efficacy and your safety.

ENG

ENG

Table of Contents

1. Introduction
2. Benefits
3. Before You Start
 - 3.1 General Description
 - 3.2 Package contents
4. Important Safety Precautions
 - 4.1 Warning
 - 4.2 Avoid to use Freequent™ if...
5. Preparing for use
 - 5.1 Charging
 - 5.2 Preparing treatment areas
 - 5.3 Patch test
6. Using Freequent™
 - 6.1 Course of treatment
 - 6.2 Before treatment
 - 6.3 Possible side effects
 - 6.4 Operating the system
 - 6.5 Post Treatment
 - 6.6 Treatment tips
7. Cleaning, maintenance and storage
8. Troubleshooting
9. Specification
10. Labeling
11. Service and Warranty
12. Copyright Reserved

1. Introduction

Congratulations on your purchase of Freequent™! Experience the radiofrequency and cooling therapy effects of Freequent™!

How does Freequent™ give you a smoother and younger skin?

Freequent™ is a new high tech home-use RF device for skin tightening. Unlike common Bi-polar RF devices in the market that only use two electrodes to deliver the RF energy to the skin, the depth of penetration of bi-polar devices for the treatment area is fixed and only just about the half the distance between the two electrodes.

However, Freequent™ uses four electrodes to deliver the RF energy to the skin. It focuses four low power beams deep into the dermis to generate a thermal effect into the skin to stimulate dermal activity by heating fibroblasts to accelerate new collagen formation and tightening collagen fibers (see Fig.1).

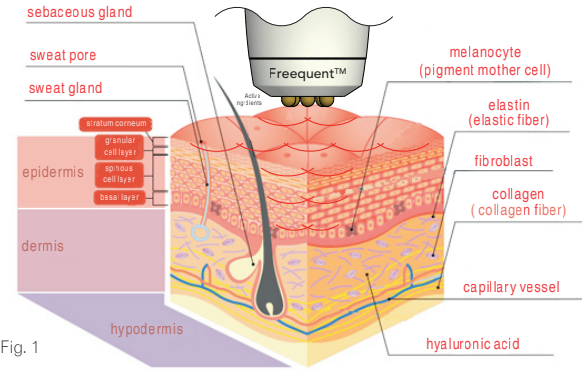
Fig. 1

As the first home-use RF device with 24 karat gold plated probes, how does this benefit your skin?

Twenty-four / 24 karat gold facial products have been used for cosmetic purposes since ancient times. This metal has been found to increase collagen production within the skin. Elasticity is improved through the restoration of elastin fibers and this dual action causes the improvement and prevention of aging in the skin.

ENG

ENG

Used as an anti-aging ingredient, it helps reduce inflammation that causes hyper pigmentation, wrinkles and age spots. It helps bring brighter-looking skin and reduce redness.

How does cool therapy benefit your skin?

The cool therapy of Freequent™ is specially designed for calming the skin after RF treatment. This professionally devised 2 min. cold treatment can reduce swelling, pores size and give a fresh feeling sensation.

2. Benefits

Freequent™ provides a professional skin care solution:

Mode	Treatment
RF for face	For cheeks, chin, the lower jaw (See Fig.2 "U1")
RF for eyes	For under and around eyes (See Fig.2 "U2")
Cool	For full face (See Fig.2 "U3")

U1-Face
U3-Cool

U2-eyes

Fig. 2

24 Karat Gold plated to maximize effectiveness

Pure 24 Karat gold plated probes provide long-term effectiveness by its unique attribute of increasing collagen production within the skin. This precious metal has been found to improve elasticity and prevention of aging in skin. (see Fig.3)

Fig. 3

High intensity energy

Freequent™ uses four probes to deliver RF energy into the skin to generate thermal effects into the dermis to stimulate dermal activity. This heating of fibroblasts accelerates new collagen formation and tightening of collagen fibers.

Designed for convenience and flexibility

With a lightweight, portable and rechargeable design, you can use Freequent™ at home or on the move.

Special probes designed for rejuvenation around the eyes

This innovative and unique tip is specially designed for delicate areas under and around eyes that need a more accurate and gentle touch. (see Fig.4)

Fig. 4

Cooling treatment

A professionally designed 2-minute cold treatment can reduce redness, tighten pores and provide a re-energised feeling of freshness.

ENG
ENG

3. Before You Start

3.1 General Description (Fig. 5)

1. Main Body
2. 24 K gold plated tip for face
3. 24 K gold plated tip for eyes
4. Button for "On/Off/ Select"
5. Light indicator:
 - White for power on
 - Green for RF mode
 - Blue for cool mode
6. Charge contact point
7. Charger Base
8. Power adaptor

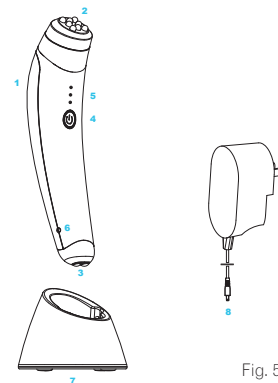
Fig. 5

Notes:

- The AC/DC adaptor is used to isolate the equipment from the supply mains.
- Charger/charger base is specified as part of the device.

3.2 Package contents

- Main device
- Charger Base
- Adaptor
- Instructions manual
- Soft Cloth
- Cleansing brush

4. Important Safety Precautions

Please read this instruction manual carefully before you use the appliance and save it for future reference. If you are unsure if Freequent™ is safe for you to use, please consult your doctor or dermatologist!

4.1 WARNING:

As with most electrical appliances, electrical parts are electrically live even when the switch is off. To reduce the risk of injury or death by electric shock:

- Do not reach for the unit if it has fallen into water or other liquid. Unplug unit immediately. Do not store charger where it can fall or be pulled into a sink or bathtub.
- Do not place or drop into water or other liquids.
- The adaptor and charger base are not waterproof. Please keep away from water.
- Keep the appliance and adaptor dry at all times.
- Do not attempt to open or repair the appliance by yourself.
- If the appliance is broken, do not touch any inner part to avoid electric shock.
- Never insert any objects into the appliance.
- Do not use the appliance if it is not working properly or if it appears damaged.
- Do not use near shortwave or microwave equipment.
- Make sure the adaptor does not get wet.
- Keep away from flammable aerosol products when used or where oxygen is administered.
- Only use parts supplied with the appliance.
- Only charge the appliance with the adaptor supplied.
- Do not use the appliance or adaptor if it is damaged.
- Do not charge when exposed to direct sunlight or heat.
- Keep this product out of reach of children. The use of this appliance by persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge can cause hazards, unless under supervision or instruction concerning the use of appliance by a person responsible for their safety.
- The adaptor contains a transformer. Do not cut the adaptor to replace with another plug as this will cause a hazardous situation.
- Turn off the appliance if it becomes unusually warm to touch.
- Switch off the appliance before cleaning.
- Stop immediately if you experience any adverse effects or develop prolonged redness.

ENG

ENG

- Freequent™ is designed for home use (as opposed to exhaustive commercial or professional use). Avoid such exhaustive use.
- Freequent™ is only designed for personal use, **DO NOT** share with others.
- Do not use any accessories, detachable parts and materials which are not described in the instructions for use.
- Do not replace the battery by yourself as it can result in an unacceptable risk.
- No modification of this equipment is allowed.
- Do not touch the output the AC/DC adaptor when using the equipment.
- Do not inhale or swallow the small parts.
- Do not interconnect this equipment with other equipment.
- This appliance is not intended for use by persons (including children) with reduced physical, This appliance is not intended for use by persons (including children) with reduced physical, given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- The battery should not replacing.

Warning

IMPORTANT! Do not use device over eyelids, lips, ears, breasts or genital areas.

IMPORTANT! Failure to follow the instructions described in this manual may cause damage or burns to your skin.

- Use this device for its intended purpose only, as described in this user manual.
- Once device is fully charged, please disconnect. Do not overcharge.
- Hands must be completely dry when operating device.
- Freequent™ RF Gel must be used during treatment according to instructions.

4.2 Avoid to use Freequent™ if...

- Under 18 years of age.
- Do not use if you are pregnant, breast feeding or nursing.
- If you have a pacemaker or internal defibrillator.
- If you have superficial metal or other implants in the desired treatment area.
- Do not use on dry or fragile skin caused by the use of chemical peels, glycolic peels, Alpha Hydroxy Acids (AHAs), Beta Hydroxy Acids (BHAs), Retin-A, Accutane and topical isotretinoin and azelaic acid.
- If you have a history of any type of cancer, or pre-malignant moles.
- If you have varicose veins in the treatment area.
- If you have pre-cancerous lesions or multiple atypical moles in the areas to be treated.
- If you have eczema, psoriasis, lesions, open wounds or active infections, such

as cold sores in the area to be treated. Wait for the affected area to heal before using Freequent™.

- If you have an impaired immune system due to immunosuppressive diseases such as AIDS and HIV, or use of immunosuppressive medications.
- If you have had surgery in the areas to be treated.
- If you have diabetes, or other systemic or metabolic diseases or congestive heart disease.
- If you have any other condition which in your physician's opinion would make it unsafe for you to be treated.
- If you have epilepsy.
- If you are under going cancer therapy.
- If you have received radiation therapy or chemotherapy treatments within the past 3 months.
- If the treatment head has cracks.
- If you have permanent fillers within the treatment area.
- If you have temporary dermal fillers within treatment area, please consult your physician before using Freequent™.
- Current use of medications, herbal preparation, vitamins or food supplements that may cause fragile or dry skin.
- Any surgical procedure, laser resurfacing or deep chemical peeling within the treatment area in past 3 months or still undergoing a process of complete healing.

Note: This list is not exhaustive. If you are not sure whether you can use the appliance, please consult a doctor.

ENG

ENG

5. Preparing for use

5.1 Charging

- Fully charge the batteries before you use Freequent™ for the first time and when the batteries are empty. Fully charging the batteries takes up to 2 hr. 30 mins
- Low battery alert: Charge Freequent™ when the light indicators are blinking. This indicates the battery is low on power.

Remarks: Due to slightly higher power consumption of the cooling feature of the device, a blinking blue light signal will indicate there is still sufficient power to use the RF function. Once a blinking green light appear, please recharge your device.

- Fully charge the appliance every 3 to 4 months, even if you do not use Freequent™ for a long time.

Directions for charging the Freequent™:

1. Switch off the appliance.
2. Insert the small plug from the adaptor into the back of the charger base and put the adaptor in the wall socket. (See Fig. 6)

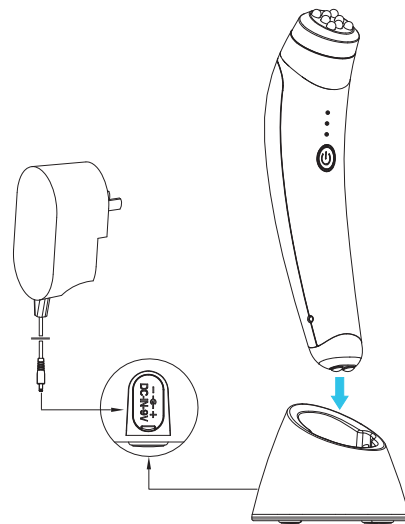
Fig. 6

3. Charging indicator light (See Fig. 7)

- The "On/Off/Mode" key illuminates in continuous red whilst charging.
- The "On/Off/Mode" key illuminates in continuous green when fully charged.

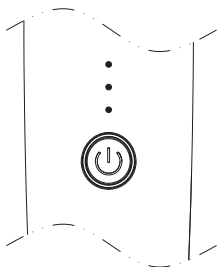
Fig. 7

Notes: It is normal when the adaptor and the appliance feel warm during charging.

Warning: NEVER cover the appliance and adaptor during charging. When charging/using, do not touch charging connector and any person simultaneously. Once device is fully charged, please disconnect power source. Do not overcharge.

4. Remove the adaptor from the wall socket and pull the small plug out of the appliance when fully charged.

5.2 Preparing treatment areas

- Start by cleaning your face with a non-oil based cleanser.
- Use the provided soft cloth to clean the surface of the attachment.

5.3 Patch test

- Familiarize yourself with the features of the Freequent™ device.
- Apply the RF gel onto the a small patch of skin, then turn on Freequent™ with the RF mode and test on a small patch of skin and wait 24 hours. Please do the same patch test with the cool mode.
- This patch test is to ensure there are no adverse reactions. If irritation occurs do not continue treatment and please consult a doctor or dermatologist.

6. Using Freequent™

6.1 Course of treatment

1. We highly recommend to follow this treatment course:
RF: 2min. for one side of face, + 2min. for the other side of face + 2mins for forehead + 2mins for under and around eyes (please use the special tip for eyes) + 2mins cool mode.
2. For optimal clinical results, use your Freequent™ device 2 times a week for a period of 6-8 weeks.
3. Please ensure you have at least a one-day interval between treatments on the same treatment area.

6.2 Before treatment

1. Cleanse the area to be treated of all creams and ointments.
2. Freequent™ RF gel must be used throughout treatment. Spread a thin, even layer over the entire treatment area.

6.3 Possible side effects

The following side effects are rare; however in the event that they do occur, they are usually temporary. However, if any of the following conditions appear and last more than 24 hours, consult your physician.

- Prolonged pain
- Damage to natural skin texture (crusting, blister, burn)
- Fragile skin
- Bruising
- Mild headache

6.4 Operating the system (Please check with section 3.1 "device description")

1. Turn on by long pressing power button. Make sure the white indicator light appears.
2. Press again for entering in RF mode. A single green indicator appears that indicates the lowest energy level.
3. Place the device on your skin, ensuring that all electronic tips are touching the skin at all times.
4. Start moving the device slowly in a circular motion, allowing the device to warm your skin to the optimal temperature. Please follow 6.1 for course of

ENG

ENG

treatment.

5. To change the energy level, press “Power” button again for higher energy level (Max. is level 3). The green light indicator will light up correspondingly.
6. The RF treatment will stop automatically after 8mins as a reminder to change to cooling treatment by pressing again “Power” button until the blue light indicators appear.
7. In order to avoid overheating the device, please power off and wait 10mins before starting the next treatment.
8. Please recharge the device if the light indicators are blinking (see Fig. 7)

Note: 10 min. is the recommended treatment time per session including 2 min. for cooling treatment.

6.5 Post Treatment

1. Turn device OFF.
2. Gently clean your skin
3. Apply a moisturizer to all treated areas.
4. Clean the device with the cleansing brush if necessary and clean with a soft, damp cloth or paper towel. Never submerge the device under water.
5. Drink more water after treatment.

6.6 Treatment Tips

1. Slow, circular and constant movements on the treatment area will allow energy to internally warm the dermis.
2. Always keep an upward movement. The treatment can be augmented by “lifting” movements. (See Fig. 8)
3. Divide the treatment area to left and right of face, forehead, under and around eyes. (See Fig. 8)

Fig. 8

ENG

ENG

7. Cleaning, maintenance and storage

⚠ Caution: Before cleaning Freequent™, make sure that the power switch is OFF and the power cord is disconnected from the charger base.

- Never clean Freequent™ or any of its parts under a tap or in the dishwasher.
- Never use abrasive cleaning agents or aggressive liquids such as petrol or acetone to clean the appliance.
- Use the soft cloth included with your Freequent™ device to clean the following parts:
 - The surfaces of the main body and charger base.
 - The surfaces of the attachments.
- Store Freequent™ in a dry place at a temperature between -10°C and 45°C.
- If you want to replace a part or purchase an additional part, please go to your local provider or contact us at sales@qlarite.com.

⚠ Caution: Protecting the environment

- At the end of its life, **DO NOT** discard Freequent™ with normal household waste. Dispose of it at an official recycling collection point.

8. Troubleshooting


Refer to this troubleshooting guide if you experience any problems with Freequent™, as this section addresses the most common problems you could encounter with this appliance and take the following precaution:

1. Ensure the electrical cord is securely connected to the charging base.
2. Ensure Power button is pressed and the light indicator appears.
3. Ensure you have applied the correct RF gel according to instructions above. Freequent™ will not work properly if the gel is not applied.
4. Ensure the probes are in contact with the skin at all times during treatment. If you still have problems, please contact your local provider or contact us at sales@qlarite.com. **Don't** open or repair the device by yourself.

9. Specification

Mode No.	QM-EF68
Lithium-ion Battery	8.4V 550mA
Charging Time	Approx. 2 hr. 30 mins
Charging Current	220mA
Adaptor spec.	Input: AC100-240V, 50-60Hz, 0.35A; Output: DC9V, 1A
Operating conditions	Temp. +15°C to +30°C
	Relative humidity: 25% to 75%
Storage and transport conditions	Temp. -10°C to +45°C
	Relative humidity: 10% to 90%

10. Labeling

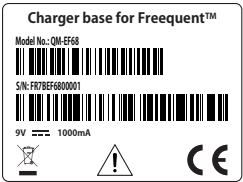ENG

ENG

Symbol	Description
	Warning/ Cautions
	Keep away from eyes
	Keep away from water
	CE Mark
	FCC Mark
	WEEE Waste Electrical and Electronic Equipment
	Refer to instruction manual
	Fragile; handle with care
	Keep the product in the dry place away from water and rain
	This way up
	Product package should be recycled

Safety Standards

- LVD: EN 60335-2-23:2003 + A1:2008 + A11:2010+A2:2015
EN 60335-1: 2012 + A11: 2014+A2:2017
EN 62233: 2008
- EMC: EN55011:2016, EN55014-2:2015, EN61000-3-2:2014,
EN61000-3-3:2013
- FCC: ANSI C63.4 & 47 CFR PART15, FCC OST/MP-5 & FCC PART 18
- Manufacturer System: ISO13485:2003 for Medical CE

ENG

ENG

11. Service and Warranty

We guarantee this product against any defects that are due to faulty material or workmanship for a 2 year period from the original date of consumer purchase. Should the product become defective within the warranty period, we will repair any such defect without charge provided there is proof of purchase. This does not mean an extension of the warranty period.

The warranty shall apply in all countries in which our product is sold via an authorized dealer.

This warranty does not cover the damage due to accidental misuse, abuse, alteration to the product or use inconsistent with the technical and /or safety instructions required.

This warranty shall not apply if the product has been dismantled or repaired by a person not authorized by us.

12. Copyright Reserved

Copyright © 2017 Qlarité SAS. All Rights Reserved.

Qlarité SAS reserves the right to make changes to its products or specifications to improve performance, reliability, or manufacturability. No part of this document may be produced or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the express written permission of Qlarité SAS. Qlarité SAS has patents and pending patent applications, trademarks, copyrights, or other intellectual property rights covering subject matter in this documents. System design and specifications may change without prior notification.

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Copyright © 2017 Qlarité SAS. Tous Droits Réservés.